

# coDiagnostiX<sup>®</sup>

## Implant Planning Software

Predictable and consistent results.

Market-proven and established solution.

300,000+ CASES


20+ YEARS



# Receive predictable results, while minimizing the number of patient visits.

coDiagnostiX is a powerful, but easy-to-use digital solution for implant planning and the design of customized and highly precise surgical guides.

- Enhanced confidence and peace of mind.
- Improved patient experience.
- Quick turnaround times.
- Reduced treatment time.
- Comprehensive case documentation.
- Advanced treatment workflows supported.

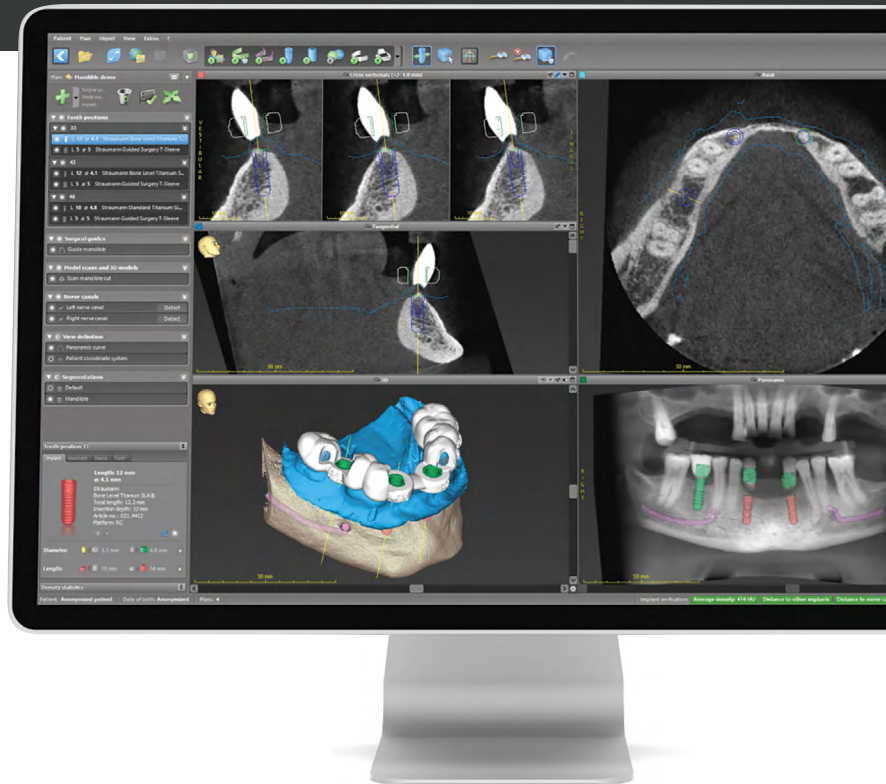

## #IChoose

“I choose coDiagnostiX because it provides a huge opportunity to improve the clinical outcomes in my clinic. It allows me to maximize performances, studying exactly where and what type of implant-prosthesis and tissue regeneration is needed.”


Dr. Andrea Pandolfi, implantologist  
Studio Dentistico Pandolfi  
Aprilia (Lazio), Italy

# A powerful and easy-to-use digital solution.

## Highly-precise and customizable surgical guides.

The coDiagnostiX guide designer offers the most flexible design tools to create surgical guides that are exactly tailored to your patient's anatomic needs.

- Time-saving guide design wizard.
- Tooth, gingiva or bone support, plus combined variants available.
- Palatal support and/or lateral fixation to facilitate edentulous cases.

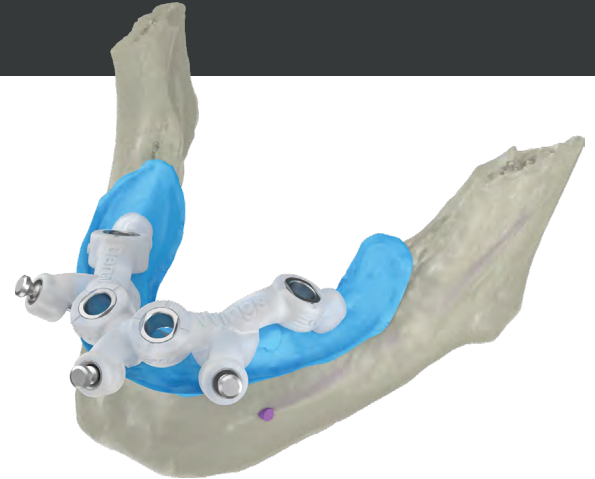

For even more flexibility, surgical guides designed with coDiagnostiX can be produced in-house, by your local lab, or through centralized service providers.

coDiagnostiX  
WORKFLOW video:



## Benefit from an open system and a completely digital workflow.

Thanks to a fully digital workflow, no scan template is required. Open interfaces ensure compatibility with all implant and sleeve systems as well as open scanners and manufacturing solutions:

- Compatible with all open (CB)CT, intraoral, model, and impression scanners, plus open CAD/CAM solutions.
- Includes a regularly updated library of implant, sleeve, and abutment systems from multiple manufacturers.
- Permits the creation of custom elements to meet individual needs.



Network with your dental partners  
efficiently using seamless  
communication technologies.

## DWOS<sup>®</sup> Connect Connecting dental professionals around the globe.

With a few clicks, scans produced with a Dental Wings 3D scanner can be sent through DWOS Connect and easily imported in coDiagnostiX.



## caseXchange<sup>™</sup> Communication platform for coDiagnostiX users.


- Send case data to a laboratory to produce surgical guide and/or provide related services.
- Exchange case data with a colleague or expert to get a second opinion or expert advice.
- Send case data to the coDiagnostiX iPad<sup>®</sup> app for case presentation to the patient.

# DWOS<sup>®</sup> Synergy

## Maximize treatment outcomes through collaborative planning.

Integrate surgical and prosthetic workflows through real-time case sharing between coDiagnostiX and Dental Wings DWOS CAD/CAM to speed up work.

For collaboration with other third-party CAD systems, a separate open interface is provided.

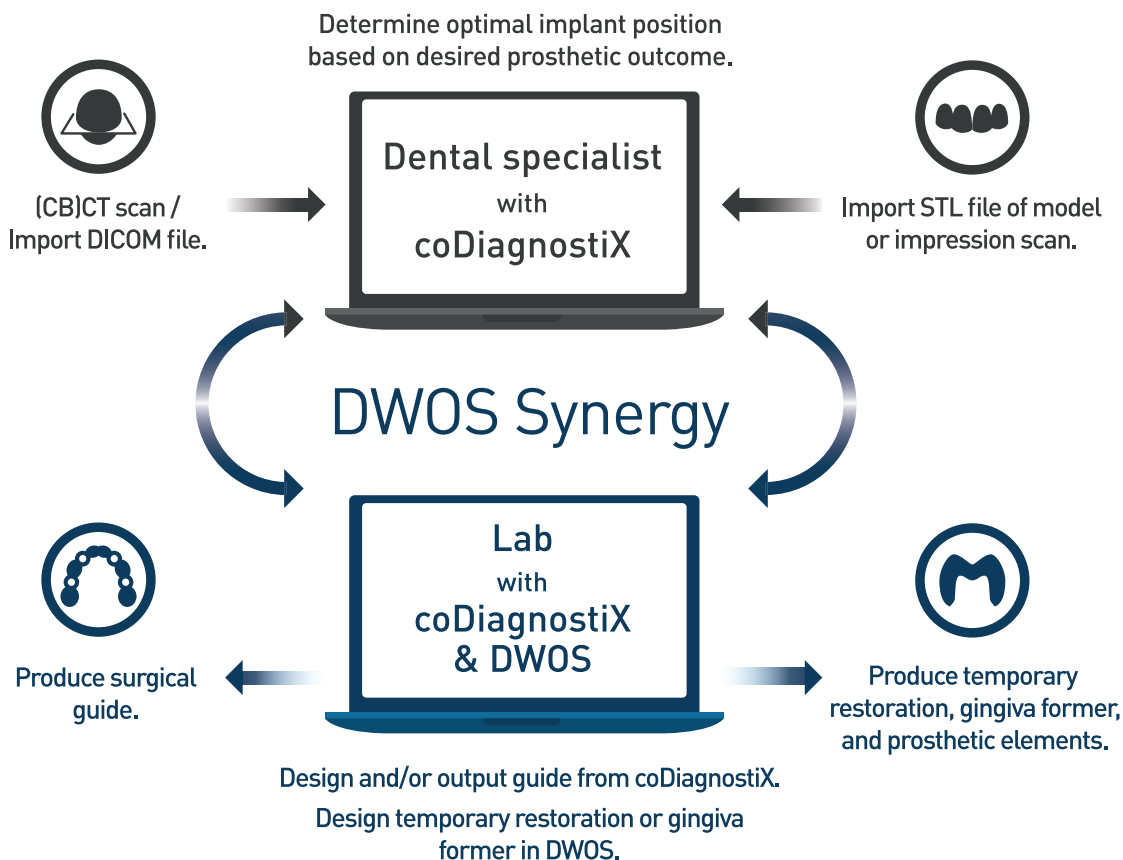



# New! coDiagnostiX EASY





Available now.

The simplified interface and full workflow from dataset alignment to export of surgical guides.

EASY mode is ideal for clinicians who outsource services to a laboratory or planning center.

- Step-by-step user guidance.
- Clear and intuitive user interface.
- Fresh and modern look.

Dental Wings offers attractive licenses with EASY-only functionality. Please refer to the [coDiagnostiX license features at a glance table](#).

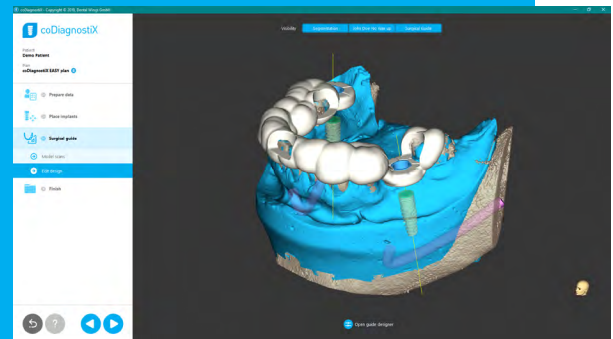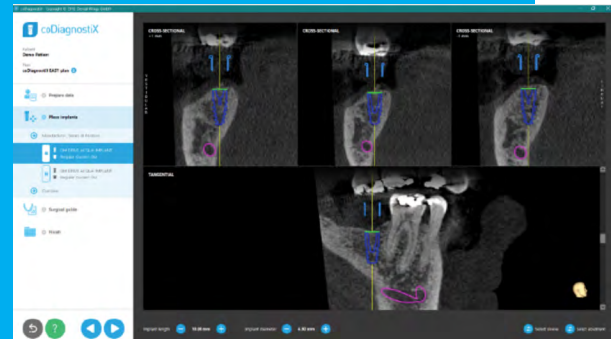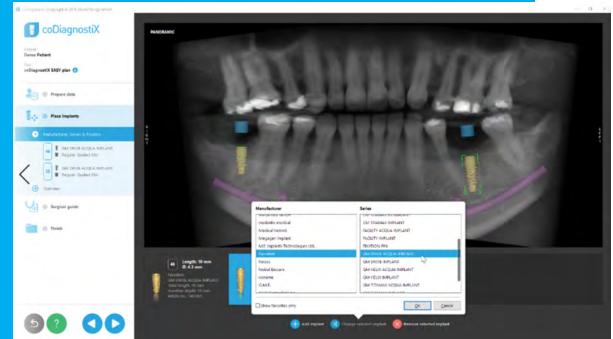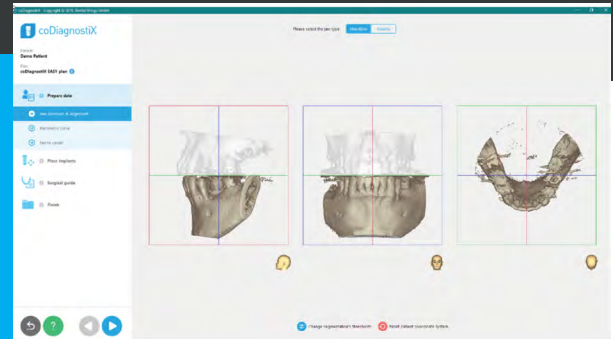

# coDiagnostiX license features at a glance

Based on coDiagnostiX 9.14

	PRODUCER	CLIENT	EASY	EASY CHAIRSIDE
<b>Planning</b>				
EASY user interface	●	●	●	●
Panoramic curve	●	●	●	●
Patient coordinate system	●	●	●	●
Planning of implants, abutments & sleeves	●	●	●	●
Comprehensive implant library	●	●	●	●
<b>Safety</b>				
Nerve detection	●	●	●	●
Distance warning	●	●	●	●
<b>Data Import</b>				
DICOM import	●	●	●	●
Model & intraoral scan import	●	●	●	●
CARES® / DWOS® Connect integration	●	●	●	●
<b>Data Output</b>				
Printouts	●	●	●	●
iChiropro	●	●	●	●
iPad viewer	●	●	●	●
External navigation system (Mininavident)	●			
<b>Surgical Guides</b>				
Surgical guide design	●	●	●	●
Surgical guide export	●			●
<b>Collaboration</b>				
Transfer cases via caseXchange™	●	●	●	
Transfer cases via archive exchange	●	●		
Multi-seat planning (network database)*	●	●		
<b>Prosthetics-driven Planning</b>				
Virtual teeth	●	●		
CARES® / DWOS® Synergy	●	●		
Third-party CAD interface	●	●		
<b>Advanced Planning</b>				
Distance measurements	●	●	●	●
Continuous distance measurements	●	●		
Angle measurements	●	●		
EXPERT user interface	●	●		
Annotations	●	●		
Augmentation	●	●		
Image management	●	●		
Implant designer	●	●		
Custom sleeve system	●	●		
Dual scan procedure	●	●		
Support for edentulous cases	●	●		
Segmentation	●	●		
Treatment evaluation*	●	●		
Virtual planning export	●			

\* Optional feature. Available upon special request.

For utmost flexibility, licenses with a limited scope of functions can be upgraded at any time.



Learn more about coDiagnostiX®: [dentalwings.com/codiagnostix-clinic](http://dentalwings.com/codiagnostix-clinic)

Download a free coDiagnostiX® demo version: [codiagnostix.com/demo](http://codiagnostix.com/demo)

Contact a sales representative:

### **Straumann North America Headquarters**

Straumann USA, LLC  
60 Minuteman Road  
Andover, MA, 01810

T 800 448.8168

F 978 747.2490

W [www.straumann.us](http://www.straumann.us)

